

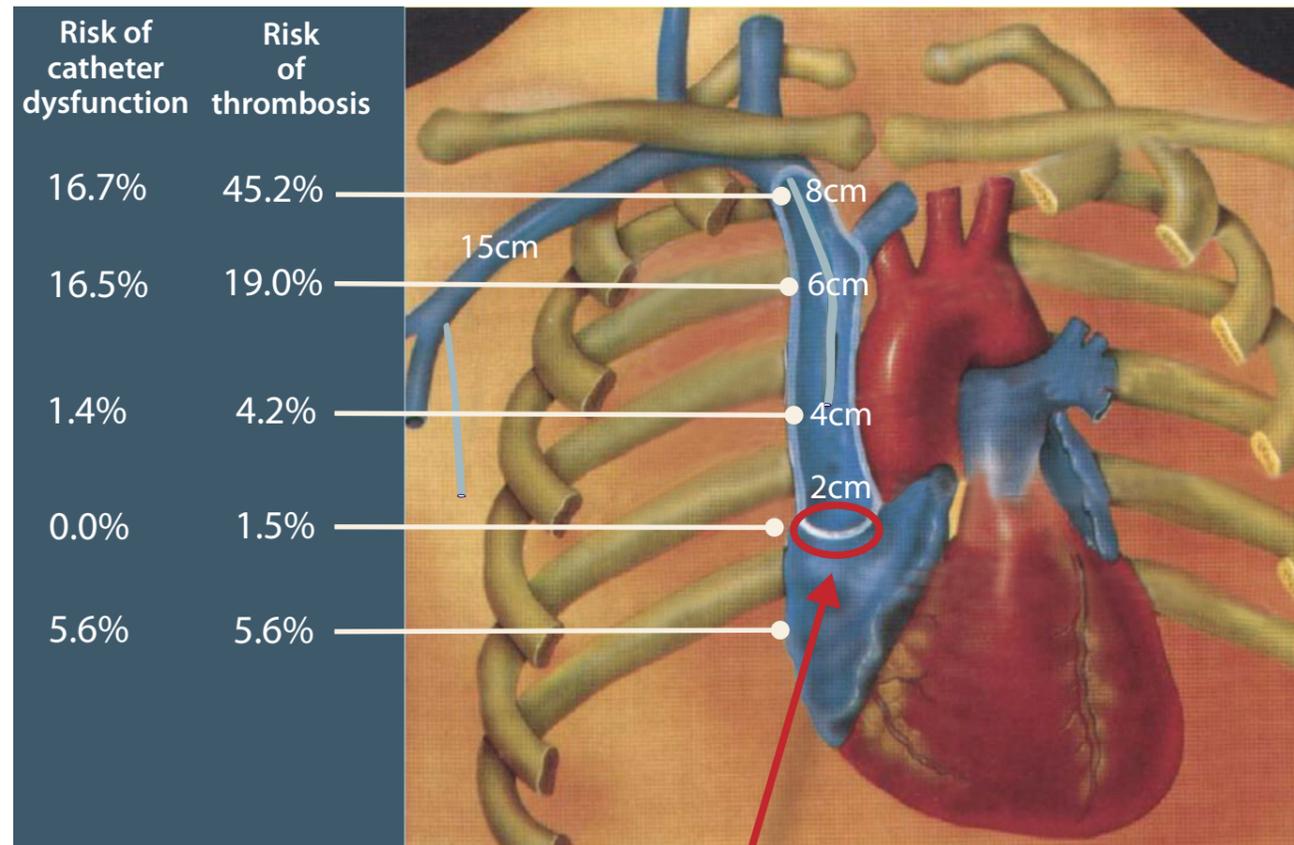
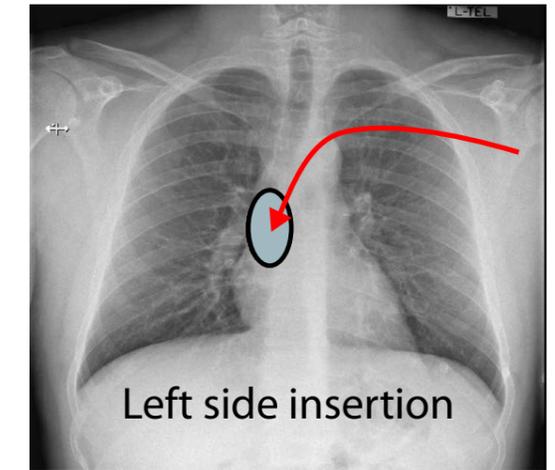


Is your PICC tip in the SAFE ZONE?

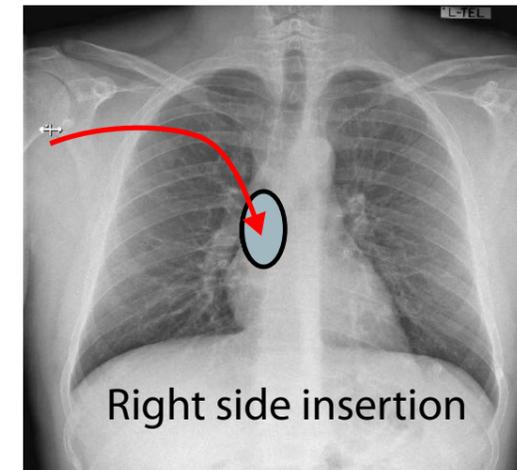

The catheter tip position with the greatest safety profile is the Cavoatrial Junction (CAJ).

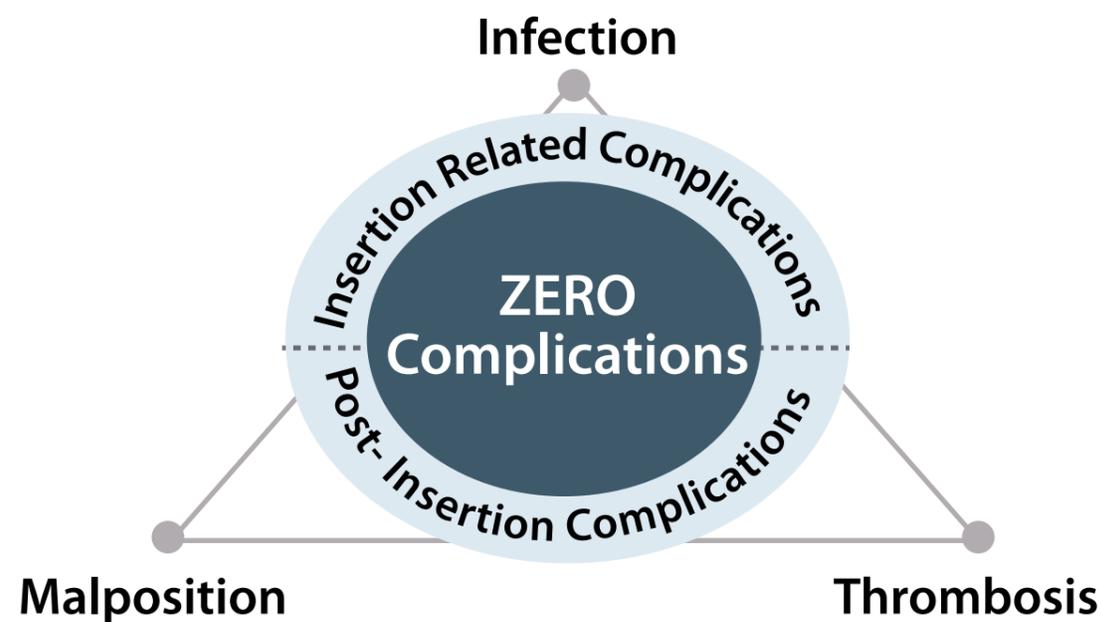

- 1** If the catheter tip *migrates out 2cms* continue to use, but monitor for any further migration.
- 2** If the catheter tip *migrates out 4 - 8cms*, only isotonic solutions should be infused. If irritant or vessicant solutions i.e. Parenteral Nutrition are required, the **PICC must be removed** to prevent complications such as vessel erosion, perforation or DVT.
Replace PICC if still required.
- 3** If the catheter tip *migrates out 8cms or greater*, **REMOVE PICC**.
Replace PICC if still required.